

POSTNATAL & NEONATAL EQUIPMENT DIRECTORY

Keeping babies warm and families close



CROYDE MEDICAL



Medical Healthcare
Furniture
In partnership with
NHS Supply Chain



The Croyde Height Adjustable Overbed Cot glides under the bed, putting the baby in front of Mum.

It allows constant and easy access, facilitating more frequent skin-to-skin contact, and promotes maternal bonding.

The mother can place the baby back in the cot without leaving the bed, reducing the risk of her falling asleep with the baby because the cot is too far away.

The Overbed Cot can be raised and lowered in one simple smooth movement to make access to the baby easier for women who have undergone caesarean sections or complicated births.

The open design makes it easy to clean, and a wide range of colours are available to suit any room theme.



HEIGHT ADJUSTABLE OVERBED COT

Product Code: BM/04/GS-KR-01ST



Chosen by birth centres, postnatal areas, NICUs, SCBUs and HDUs, the Croyde Height Adjustable + Storage cot offers the convenience of customisable and practical under-cot storage.

The large cupboard is accessible from both sides of the cot with double soft-close doors, ensuring it is suitable for small rooms or tight spaces.

Height adjustment enables the cot to be raised or lowered to the desired height to limit bending or stooping, reducing musculoskeletal strain on the midwife's back and shoulders.

You can specify an optional foldable shelf or oxygen cylinder holder, and a wide range of colours are available to suit any room theme.



HEIGHT ADJUSTABLE STORAGE COT

Product Code: BM/05/GS



AQUATHERM WARMING MATTRESS

Product code: 7674-11282A73

The Croyde Aquatherm Warming Mattress is designed to provide controlled heat to newborn infants, particularly those who are premature or have difficulty maintaining their body temperature.

The mattress uses a fire retardant foam/ rubber and 'Aquagel' technology to deliver a non-water based warming system. A heating element is incorporated for safe and evenly spread heat as well as two highly-sensitive sensors for fail-safe temperature control.

When disconnected, the Aquagel technology retains residual heat to create a stable warming environment during transport.

The control unit has user-friendly settings to regulate the heat output of the mattress and also contains an alarm system for enhanced patient safety.



MEDWARM NEONATAL WARMING MATTRESS

Product Code: IM-MS

The Medwarm Neonatal Warming Mattress sits in infant cots, offering controlled heat to keep newborns warm and comfortable - particularly premature babies and those that struggle to maintain their body temperature.

Medwarm stays fully operational when unplugged thanks to its rechargeable battery. The robust control unit will operate for up to four hours, keeping the mattress actively heated and reducing the risk of hypothermia.

The mattress operates at a temperature of 39°C, reaching 37°C in just 5-7 minutes, helping to protect an infant's core body temperature and avoid hypothermia after birth. The temperature is easily controlled, with an accuracy of 0.1°C.

The mattress core is made using memory foam, ensuring maximum comfort for baby whilst remaining lightweight for easy management and transport.



MAKE THIS A HOT COT

Add a Croyde Cot to create a fit for purpose Hot Cot. The Croyde Hot Cot provides a height adjustable, ergonomic and fit for purpose solution for neonatal warming.

Croyde Overbed Hot Cot with Aquatherm CAP-CHC01

Croyde Storage Hot Cot with Aquatherm CAP-CHCws-CSTD



MAKE THIS A HOT COT

Add a Croyde Cot to create a Hot Cot. The Croyde Hot Cot provides a height adjustable, ergonomic and fit for purpose solution for neonatal warming.

Croyde Overbed Hot Cot with Medwarm Neo - CAP-CHCM01 - M55

Croyde Storage Hot Cot with Medwarm Neo - CAP-CHCM03 - M55





Perfect for birth centres and postnatal & neonatal areas, the Croyde Next To Me cot is height adjustable and features a lowering side, bringing mother and baby closer.

The cot features intuitive release buttons on both sides, making it quick and easy to lower each side panel for hassle-free access. Engineered with your baby's safety in mind, both release catches must be pressed simultaneously, ensuring secure use at all times.

When lowered, the side panels offer full-length access to the cot, making it ideal for bedside positioning. This design enables mothers to comfortably and safely reach their baby for night-time feeds and bonding, without the strain of lifting. It's especially helpful for mothers recovering from a C-section, as it eliminates the need to bend and lift their baby from the cot.



CROYDE NEXT TO ME COT

Product Code: BM4-1



The Croyde Twin Cot features a spring-assisted height adjustment pedal, allowing for effortless raising and lowering in one smooth motion. This allows clinicians to adjust the height for neonatal examinations and mothers easy access to their newborns.

The easy-to-use release buttons on each side allow the side panels to be lowered for convenient access.

Designed with safety in mind, both catches must be released before lowering each side panel.

When lowered, the side panels provide full access along the side of the cot and enable it to be positioned next to the bed, allowing mothers to safely reach their babies.



CROYDE TWIN COT

Product Code: SBS-2.0-Twin-12



LIFETHERM 2000 WARMING BED

Product Code: CAP-LFTH2000ha

The Lifetherm 2000 provides a safe, supportive environment for infants up to 8kg, easing the transition from incubator care. From first breath to full stability - Lifetherm 2000 ensures every newborn starts warm, safe and strong.

The Lifetherm 2000 supports early stabilisation by maintaining precise and reliable thermoregulation, helping to minimise avoidable transfers to SCBU or NICU.

As well as even heat distribution, heat retention can be enhanced with an optional lid - creating a controlled environment that reduces escalation of care while preserving intensive care resources.

Access is vital when caring for fragile newborns. Designed to give both parents and clinical teams immediate, unhindered access to the baby, all four sides fold down for clear, barrier-free reach.

Both set and actual temperatures are clearly displayed, and instant alerts to any changes helps teams respond quickly and confidently in even the most critical situations.



LIFETHERM 2003 WARMING BED

Product Code: CAP-LFTH2003ha

The Lifetherm 2003 Warming Bed combines neonatal warming with critical care support, featuring an overhead radiant warmer. It ensures safe, reliable thermoregulation while streamlining workflows and improving outcomes for vulnerable infants.

The Lifetherm 2003 Warming Bed supports early stabilisation by delivering precise, reliable thermoregulation, helping to reduce avoidable transfers to SCBU or NICU.

The integrated overhead radiant warmer provides immediate thermal support during critical interventions. Optional taller side panels and a lid are available to further improve heat retention.



Standard Specification - Lifetherm 2000

- A choice of low or high side panels
- Height adjustment via foot pedals
- Heated surface (28-38.5°C, adjustable in 0.1°C steps with automatic safety controls)
- Silicone gel mattress (choose with or without safety panes)
- Audio-visual alarms

Standard Specification - Lifetherm 2003

- All items listed above plus:
- Integrated overhead radiant warmer
- Side poles

Optional extras

- Medical rail on three sides - 0215570
- Drawer block (includes attachment set) - 0125510
- Folding plexiglass hood (only to be used with high side panels) - 0215550
- Non-folding plexiglass hood (only to be used with high side panels) - 0215560
- Skin temperature sensor (only available for the Lifetherm 2003) - 0216021



ISOTHERM RADIANT WARMER


Product Code: 0210080hul100

The Isotherm Radiant Warmer helps to create a warm environment and maintain comfort. With adaptable installation, user-friendly heat control and a clinician-focused design, it provides a reliable radiant warming solution across perinatal care settings.

Four Clearly defined heat settings make temperature adjustment simple and intuitive, reducing confusion and ensuring consistent care across staff.

Developed with clinicians in mind, Isotherm combines practical functionality with patient safety. Its integrated examination light offers both dim and bright settings for clear visibility without disturbing the infant, while the user-friendly interface allows reactive adjustments in busy clinical settings.

Simple and easy to use control panel

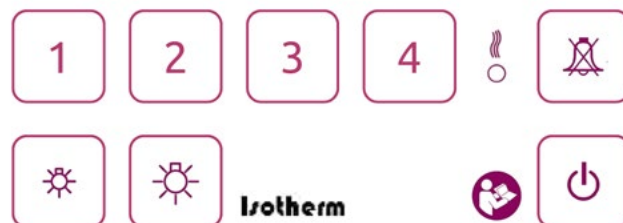

BILLIBLUE LIGHT - PHOTOTHERAPY

Product Code: 0218020

The Bilibluelight Phototherapy Lamp is a compact, lightweight solution for the effective treatment of neonatal jaundice. Combining advanced LED technology with a user-friendly design, it helps healthcare professionals deliver safe, reliable care in busy neonatal settings.

Designed to support clinical teams with fast, reliable operation, a single button control ensures consistent phototherapy delivery, minimising training needs and eliminating staff confusion.

Ultimately, its lightweight, modern design makes it easy to handle, reposition and store - helping staff save time and work more efficiently during neonatal jaundice treatment.


ELIOS

Product Code: 8037

ELIOS is a multi-functional system offering 3 functions from one mobile light head. It delivers intensive phototherapy with overhead radiant warmth and white light examination.

ELIOS is a comprehensive solution for neonatal care and eliminates the need for multiple, mobile lights. The lamp is designed for safe and easy placement over cribs and changing tables.

ELIOS has a robust, height-adjustable steel stand and an anti-tipping base equipped with swivelling wheels for smooth manoeuvrability.

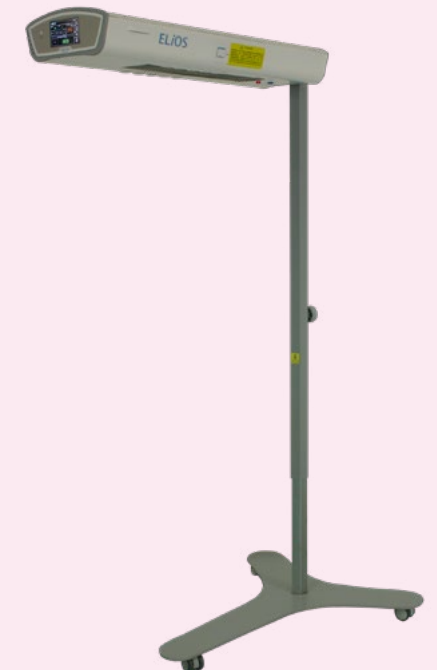



Available in a standard or compact size, the Croyde Height Adjustable Baby Assessment Table allows clinicians to position the table at the perfect height for neonatal examinations and baby changing.

Depressing the foot pedal enables staff or parents to adjust the table to the ideal position and protect them from physical stress and strain. NIPE screening becomes instantly more comfortable and stress-free.

The high-height position makes cleaning easier, while the low-height position makes the table accessible for all heights of carers, including wheelchair users. It is lightweight and easy to move between bed spaces, with brakes to secure its location.

Choose from a storage cupboard with shelving or a drawer unit. Both units have soft-close doors and ample space for all items.



The Croyde Baby Bath is a height adjustable baby bath with out Easy-Pour Tub for quick and simple emptying.

The tub can be raised and lowered using the foot pedal to tailor the height to suit any care giver and improve ergonomics when using or emptying.

The Easy Pour Tub design features a pointed spout to make emptying the bath efficient whilst reducing spillages. The handle is easy to grip when tipping and can also be used as a small shelf to store soap sachets, cotton wool etc when in use.



HEIGHT ADJUSTABLE BABY ASSESSMENT TABLE

Standard Product Code: SPTM-02-D-HA
Compact Product Code: SPTM-02COMPACT-D-HA

CROYDE BABY BATH

Product Code: CAP-TBB01



CROYDE KANGAROO RECLINER CHAIR

Product Code: TX15

The Kangaroo Recliner is an ideal choice for Maternity and Neonatal settings, combining full-recline comfort with smooth patient transfer functionality and a versatile selection of optional accessories to support both care and recovery.

The Croyde Kangaroo Recliner stands out with its unique ability to recline completely flat - a game-changer for exhausted parents in need of real, restorative rest.



Unlike many traditional recliners that can only partially recline - which can place unwanted pressure on sensitive areas, such as perineal or C-section wounds - this fully flat position provides true comfort and relief. A fully reclined position is also ideal for caregivers to perform antenatal palpation and take symphysis fundal height (SFH) measurements.

Its ergonomic push handle allows caregivers to move patients smoothly and securely between clinical areas in an emergency, reducing strain and improving workflow.



CROYDE MATERNITY & NEONATAL RECLINER

Product Code: ASTM

A multi-use bedside package, perfect for maternity wards. The recliner facilitates ideal infant feeding positions for mother and baby. The reclining backrest allows babies to feed in the semi-recumbent chest-to-chest position that is most comfortable for the mother. Padded armrests give comfortable arm placement and support the chosen feeding position. For partners and carers, the recliner positions them perfectly for skin-to-skin contact when bonding with the baby.

The padded foot section of the Croyde Maternity and Neonatal Recliner gently raises the level of the feet closer to the heart, improving lymphatic drainage, and reducing swelling in the feet, ankles, and legs after birth. By elevating the feet, the lower back also benefits from reduced strain.

The chair also allows overnight stays for partners without needing additional equipment, replicating the home environment and improving bonding with birth partners.

The Croyde Infant Feeding Pillow and Footstool are also available to purchase separately.

Available in a range of colours.





CROYDE ELECTRIC MATERNITY RECLINER

Product Code: GeoSTIS20/GeoSTIS25

The Croyde Electric Maternity Chair aids sitting, reclining, and rising at the touch of a button, and supports women who have undergone invasive surgeries. The reclining backrest allows babies to feed in the semi-recumbent chest-to-chest position that is most comfortable for the mother. Padded armrests give comfortable arm placement and support the chosen feeding position.

It makes the immediate postnatal period more comfortable for women who have restricted mobility, who are experiencing pelvic girdle pain or are recovering from invasive procedures.

The Croyde Electric Maternity Recliner Chair facilitates overnight stays for partners without additional equipment needing to be brought into the ward. This helps to replicate the home environment and improve bonding with birth partners.



CROYDE SLEEPER CHAIR

Product Code: STARGAZER

Perfect for use in hospital settings where overnight stays for partners/ dependants are required. The Croyde Sleeper Chair is easily deployed from chair to comfortable, full-length bed in just a few seconds. It also features a chaise longue position. Integrated wheels make moving simple and lightweight and it is designed to be wiped-down between uses in the hospital setting.

Available in a range of colours.



"The sleeper chairs are both spacious and comfortable, whilst not taking up too much room in our side rooms. They have helped us to ensure that our patient's relatives have a better experience during overnight stays, which are often in very difficult circumstances"

Inpatient Matron at Basildon & Thurrock, University Hospitals NHS Foundation Trust






Our electric profiling bed features hi-lo, backrest and thigh section adjustment via a back-lit handset control. It is available in 120cm or 140cm widths, or as a 90cm single bed if space is tight.

The bed has a weight capacity of 315kg and has integrated CPR release for emergency situations.

It is supplied with a hospital-specification mattress and can be cleaned with approved wipes and disinfectants. Choose from a range of lacquered wood finishes to create a home-from-home environment.

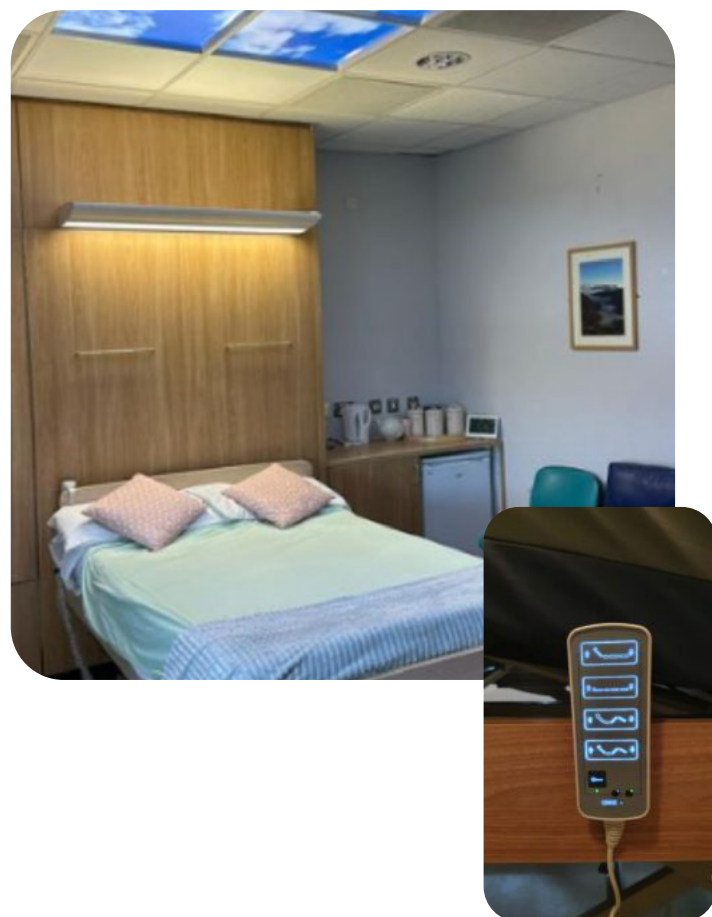

The Croyde Wall Bed is available as a single bed or a double bed and can be easily raised and lowered with just one hand.

The beds can be purchased independently or with alongside cabinetry such as wardrobes, cupboards, shelves or storage boxes.

Integrated lighting, USB charging points and power options are optional and a comprehensive delivery and installation service is available. The cabinetry can be adapted to accommodate gas systems, nurse call systems or other fixed wall systems.

A wide range of colours and finishes is available and each bed is supplied with a hospital-specification mattress that fits perfectly to the frame and easily folds away with the bed.

All items can be cleaned with hospital-approved wipes and chlorine-based solutions.



CROYDE BEREAVEMENT SUITE BED

Product Code: XB1P8ABD106

CROYDE WALL BED

Product Codes:

Double wall bed: CDWB

Double wall bed & cabinetry: CDWC+C

Single wall bed: CSWB



MIMMAMÀ INFANT FEEDING CHAIR

Product Code: 365900 S/N



The MimmaMà is specifically designed for infant feeding to accommodate a variety of different positions including cradle, cross and rugby ball options.

The shape of the chair helps to maintain the correct position of the child on the nipple, avoiding fissures. It also allows the mother to alternate feeding positions, freeing the milk ducts and avoiding mastitis.

Ideal for twin feeding due to the dual backrest spaces. There is excellent support for the lower back, arms and shoulders.

The Croyde Infant Feeding Footrest and the Croyde Infant Feeding Pillow can be used alongside the chair to provide further support for mother and baby.

Available in a range of colours.



CROYDE FEATURE WALL PAPER

Product Code: Feature Wall (Wallpaper)

Bring colour to any hospital space using the Croyde Feature Hospital Wallpaper. A custom made wallpaper is an ideal way to cover a large area and make a great impact.

There is a wide range of designs available, from photographs of a local area to a calming coastal scene or landscape. These uplifting designs can transform clinical areas into welcoming spaces and help to create oxytocin-boosting birth environments.

Usually taking just a couple of hours to complete, this feature wallpaper can be installed without causing too much disruption to normal hospital routines and appointments.



BREAST MILK FRIDGES

Product Code: BWR125

The Croyde Breast Milk Fridges are available in a range of sizes, from sixty litre capacity to four hundred litre capacity.

The digital control panel allows for the temperature selection, between 1-5 °C, with a single press 3 °C, setting for the optimum breast milk temperature. Audible and visual temperature alarms immediately alert healthcare practitioners of any major temperature fluctuations and help to ensure the quality of the milk.






QED BLANKET WARMING CABINETS

Product Code: BW-ECO

The QED Blanket Warming Cabinets are available in a range of sizes accommodating 4-50 cotton blankets in baskets. These cabinets can be used in range of settings, including birth environments, postnatal & neonatal, operating theatre recovery, critical care, A&E and other hospital wards where rapid patient warmth is important.

All QED Warming Cabinets are equipped with PID digital temperature controls. Blanket warmers have a maximum temperature of 80°C but can be set to a desired temperature.

The chamber is made from medical grade stainless steel and the external of the cabinet is powder painted to allow for easy cleaning.

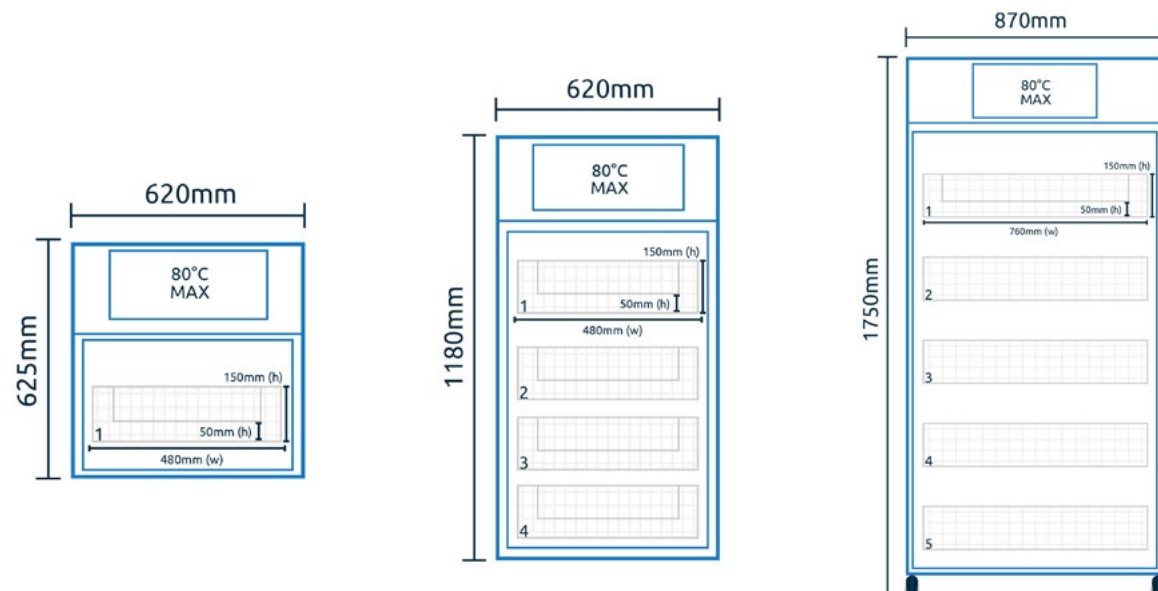

QED CONTROLLED ROOM TEMPERATURE CABINETS

Product Code: CRT

Offering specialised storage designed to maintain ambient room temperatures, the QED CRT Cabinets ensure safe storage of medications for use in clinical and pharmaceutical settings.

For medicines required to be kept under 25°C, ensuring room conditions are kept at an ambient temperature is essential for maintaining their efficacy. If temperatures get too low or high, the potency and shelf-life of these medicines can reduce significantly.

As well as storing such medicines, the QED CRT Cabinet keeps them at an optimal temperature. Not only does this ensure their efficacy, but it maintains their shelf-life - saving them from going to waste.



QED CONTROLLED DRUGS CRT CABINET

Product Code: CD-CRT

QED Controlled Drug CRT cabinets provide top security and compliance for storing controlled substances. They also feature integrated temperature control to maintain drug efficacy and reduce waste.



Solutions Beyond the Incubator

Keeping babies warm and families close.

Continuing Care



Overbed
Cots



Storage
Cots



Advanced
Care


Practical Warmth



Underbody
Warming



Radiant
Warmth

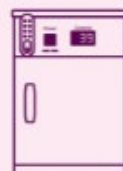

Warming
Cabinets

Comfort & Connection



Recliner
Chair



Sleeper
Chair



Kangaroo
Recliner



CROYDE MEDICAL



Scan to discover
the range