

MEDWARM PATIENT WARMING SYSTEM

Offering controlled heat to help patients maintain comfortable body temperatures in theatre rooms, reducing the risk of hypothermia.



The Medwarm Patient Warming System can be made up of the Medwarm Mattresses, Medwarm Blankets and Medwarm Control Units in your chosen configurations.


It offers efficient and intuitive conductive heating, ensuring patients' body temperatures are effectively regulated.

The mattress operates a temperature range of 30-40°C, which is adjustable in 0.1°C increments for precise control and accuracy. Heat is then distributed evenly and consistently throughout the mattress.



SIMPLE & DURABLE DESIGN

The Medwarm W-300 Control Unit offers a perfect blend of simplicity and durability. With a user-friendly interface, the W-300 features a modern, robust design built to withstand frequent use and the demanding conditions of hospitals.

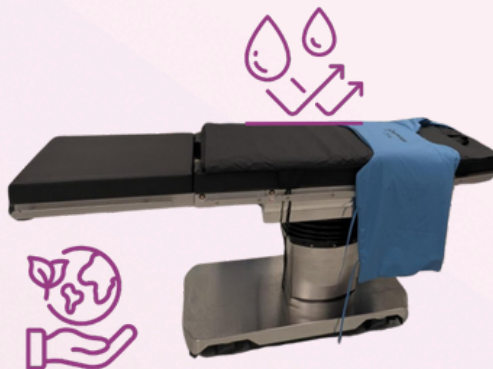

SUSTAINABLE SOLUTION

The Medwarm Mattress and Blanket provide a reusable and eco-friendly solution for patient warming that doesn't use consumables. With a 10-year lifespan, they offer extensive usage before replacement is required - offering unmatched longevity.



INFECTION CONTROL

The Medwarm Mattress and Blanket are perfect for environments where infection control is critical. They use conductive heat - eliminating air movement and minimising the risk of disturbing sterile fields or spreading airborne contaminants.



COMFORT & SAFETY PRIORITISED

Built-in alarms make up a comprehensive security system that detects and alerts if the mattress or blanket gets too hot, malfunctions, disconnects or short-circuits.

The mattress uses X-ray-permeable carbon fibre, allowing patient warming during scans, while memory foam adds comfort and helps prevent bed sores during long procedures.

EVERYTHING YOU NEED TO KNOW

Click the links below to access
the Croyde Knowledge Hub

[How to Buy](#)

[Dimensions & Specifications](#)

[Delivery & Implementation](#)

[Warranty & Servicing](#)

[Cleaning Instructions](#)



CROYDE MEDICAL

CONTACT US

Tel: 01367 241194

Email: info@croydemedical.co.uk

[**www.croydemedical.co.uk**](http://www.croydemedical.co.uk)

Croyde Medical, Unit A4 GBP1, Grove Business Park,
Downsview Road, Wantage, OX12 9FA